

Enteric Range CE

Detection of

Rotavirus

Adenoviruses (all groups)

Adenoviruses (group F, serotypes 40 & 41)

Cryptosporidium

Giardia

***Escherichia coli* O157**

Rotavirus and Adenovirus (all groups)

Rotavirus and Adenoviruses (group F, serotypes 40 & 41)

Cryptosporidium* and *Giardia

Respiratory Range CE

Detection of

Respiratory Syncytial Virus

Respiratory Adenovirus

Influenza A Virus

Influenza A&B Virus

Biomolecular Range CE

Detection of

***Leishmania* spp.**

Trypanosoma cruzi



Rapid antigen diagnostic tests

Rapid Nucleic Acids Amplification Diagnostics

R & D « on demand »

For more scientific
or technical information,
visit our website:

<http://www.corisbio.com>

Innovative solutions for more effective diagnostics



A quick and reliable answer
for diagnostic labs, emergency centres
and doctors' offices



Science Park CREALYS

Rue Jean Sonet 4A - 5032 Gembloux - Belgium

Ph.: +32 (0)81 719 917 - Fax: +32 (0)81 719 919

info@corisbio.com

Distributed by:



Company Profile

Coris BioConcept is a high technology company in the biomedical field.

The main activity of Coris BioConcept is the production and marketing of diagnostic kits for the detection of human enteric and respiratory pathogens.

Tests are performed on coated nitrocellulose strips and results are observed with colloidal gold particles or latex microspheres.

Once the strip is dipped into diluted sample, liquid samples and conjugate(s) migrate by capillary movement. At the end of the migration process, one or two coloured lines appear indicating either a negative or a positive result respectively.


Coris BioConcept works in a QC/QA frame.

The quality policy has been certified ISO 9001:2000 and ISO 13485:2003 by SGS.

Moreover all Coris BioConcept's products fulfil the CE marking requirements.

Advantages of Coris BioConcept rapid antigen diagnostic tests:

- Accurate
- Specific and sensitive
- One-step test
- Minimal sample manipulation
- Simple and easy-to-use
- No special training needed
- Results in less than 15 minutes
- No special equipment required

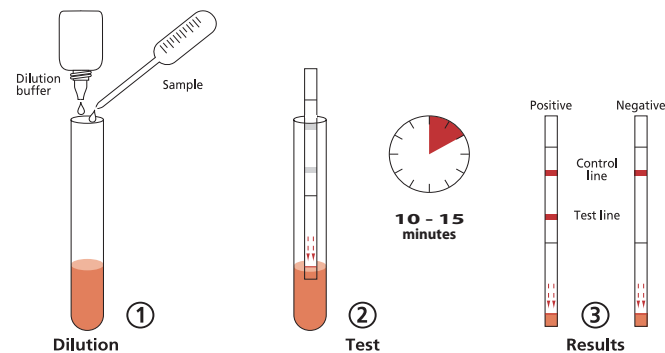

Product range

Rota-Strip
Adeno-Strip
40/41 Adeno-Strip
Combi-Strip
Crypto-Strip
O157 Coli-Strip
Gastrovir-Strip

Giardia-Strip
Crypto/Giardia Duo-Strip
RSV Respi-Strip
Adeno Respi-Strip
Influ-A Respi-Strip
Influ-A&B Respi-Strip

Strip Range

Packaging of 25 tests (or 20 tests)

- Each kit contains:
 - One tube with 25 or 20 strips
 - Dilution / extraction buffer
 - Instructions for use



Uni-Strip Range

Packaging of 10 tests

- Each kit contains:
 - Individually packed tests (pouches)
 - Assay tubes
 - Dilution / extraction buffer
 - Sampling material
 - Tube rack
 - Instructions for use

