



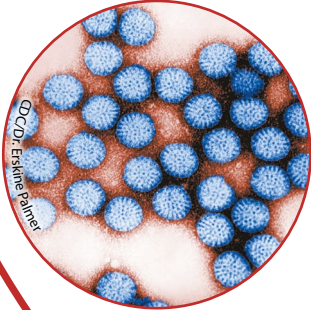

Combi-Strip


Rapid diagnostic test for *in vitro* detection of Rotavirus & Adenovirus in stool specimens

- Easy to use
- Easy to read
- Specific
- Sensitive
- Reliable
- Rapid (10')

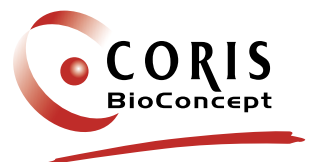

www.corisbio.com

Combi-Strip

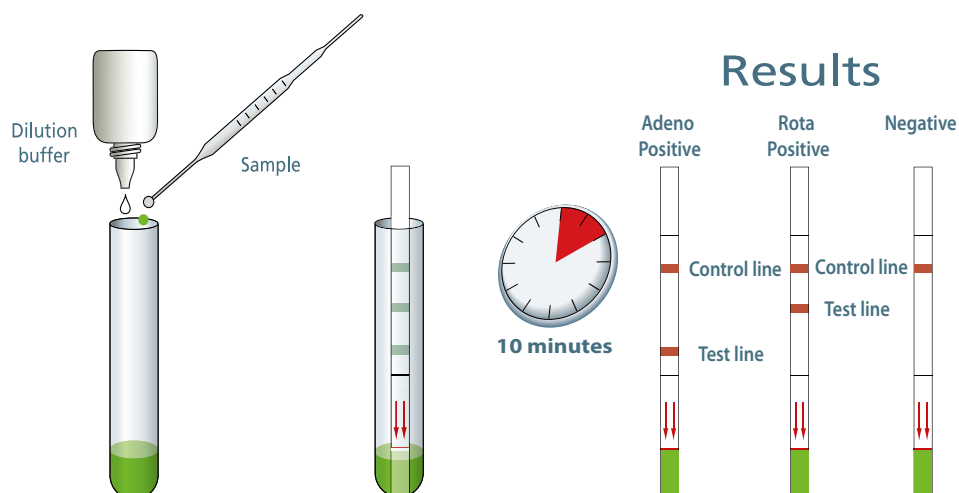

Performance of the kit


	Compared to ELISA Rotavirus detection (n=214) (95% confidence interval)		Compared to ELISA Adenovirus detection (n=130) (95% confidence interval)	
Sensitivity	99.1 %	(94.1 to 100 %)	100 %	(69.9 to 100 %)
Specificity	100 %	(95.7 to 100 %)	100 %	(96.1 to 100 %)
PPV	100 %	(95.6 to 100 %)	100 %	(69.9 to 100 %)
NPV	99.1 %	(94.3 to 100 %)	100 %	(96.1 to 100 %)
Agreement	99.5 %	(213/214)	100 %	(130/130)

For more technical information

client.care@corisbio.com

Ordering Information

PRODUCT NAME	DESCRIPTION	STORAGE	CODE
Combi-Strip	25 tests, buffer	4°C to 30°C	C-1004
Combi K-SeT	20 cassettes, buffer	4°C to 30°C	K-1204
Combi K-SeT with collection set	20 cassettes, 20 faecal sampling systems	4°C to 30°C	K-1504

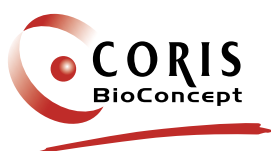

CREALYS Science Park
Rue Guillaume Fouquet 11
5032 Gembloux - BELGIUM
T +32 (0)81 719 910 - F +32 (0)81 719 919
sales@corisbio.com - www.corisbio.com

— Distributed by —



MADE IN BELGIUM

D-CUS-10-04-5/09
March 2023